



Obturation Success for Everyone



Introducing BioRoot™ Flow Bioceramic Root Canal Sealer



Easy

Ready-to-use syringe provides easy extrusion directly into canals and you maintain your preferred obturation technique.



Successful

Highly-pure calcium silicate mineral formulation (resin-free) means no shrinkage⁽¹⁾, limited bacterial growth, and increased biocompatibility and bioactivity which reduce the risk of adverse reactions.



For Everyone

Designed for cold and warm techniques with excellent flow into accessory canals⁽²⁾, tight adhesion to dentin, and gutta-percha, and reproducible results in all root canal treatments from simple to complex.⁽³⁾

Septodont

205 Granite Run Dr., Suite 150
Lancaster, PA 17601

800.872.8305

www.septodontusa.com

(1) Internal data. Pr Colon report

(2) Drukteinis & Camilleri : Bioceramic materials in clinical endodontics. Springer 2020


(3) BioRoot™ flow SEM interfaces. Dr Camilleri



Technical Insights

Innovative tip allows direct and precise placement

- ▶ Flexible tips allow bending and access to all root anatomies
- ▶ Precise 21-gauge diameter for optimal delivery limiting the risk of going too deep

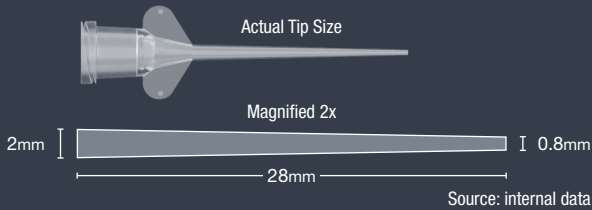


Highly Radiopaque

- ▶ 6.3mm Al radiopacity
- ▶ Easily visible on X-rays


Easily Retreatable

- ▶ Easy and fast to remove if necessary
- ▶ Use hand files or rotary NiTi Files
- ▶ BioRoot Flow is removed to working length in less than 10 minutes.

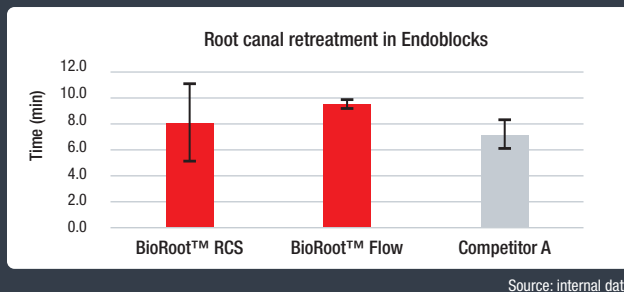

Long Lasting pH

- ▶ High pH is maintained over time creating an alkaline environment in which bacteria will not grow

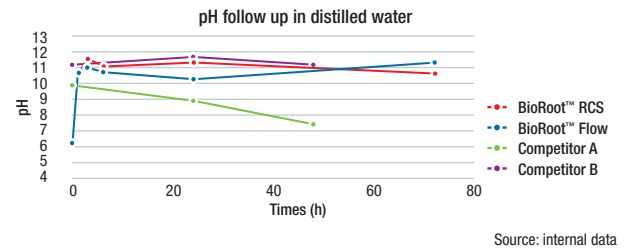

Technical Features

- ▶ Working time | > 60 min
- ▶ Setting time | 5 to 6 hr

Source: internal data; Dr Camilleri

Product Information

- ▶ 1 x 2 syringe
- ▶ 20 intra-oral tips
- ▶ 1 x finger grip
- ▶ Item #01E0510



Septodont

205 Granite Run Dr., Suite 150
Lancaster, PA 17601

800.872.8305
www.septodontusa.com



R0322-1K