

NEW!

The Light That Moves With You



Demi™ Pro Curing Light

Kerr™



Cure With Control, Not Fatigue

- ✓ Ergonomic comfort
- ✓ Pen-like precision in a lightweight, balanced body at just 115g
- ✓ 360° access with a low-profile tip



"The Demi Pro curing light is very lightweight and made curing at every angle comfortable."

- Dr. Shahira Saad



The Demi Pro Advantage



All-Day Comfort

Reduce fatigue with the lightweight, ergonomic design.



Easy Access

Reaches every angle with the unique 360° and low profile rotating tip.



Deeper Cures

Delivers consistent curing performance at clinical distances, thanks to its collimated output and improved beam profile.



Safety

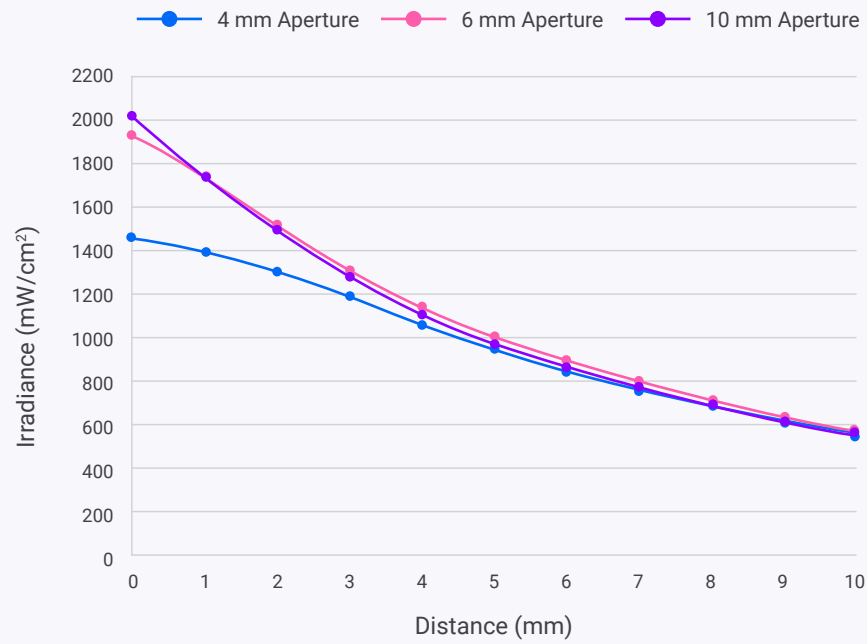
Minimize the risk of overheating with patented PLS (Periodic Level Shift) Technology and an integrated radiometer.



360°
PRECISION
REACH EVERY ANGLE

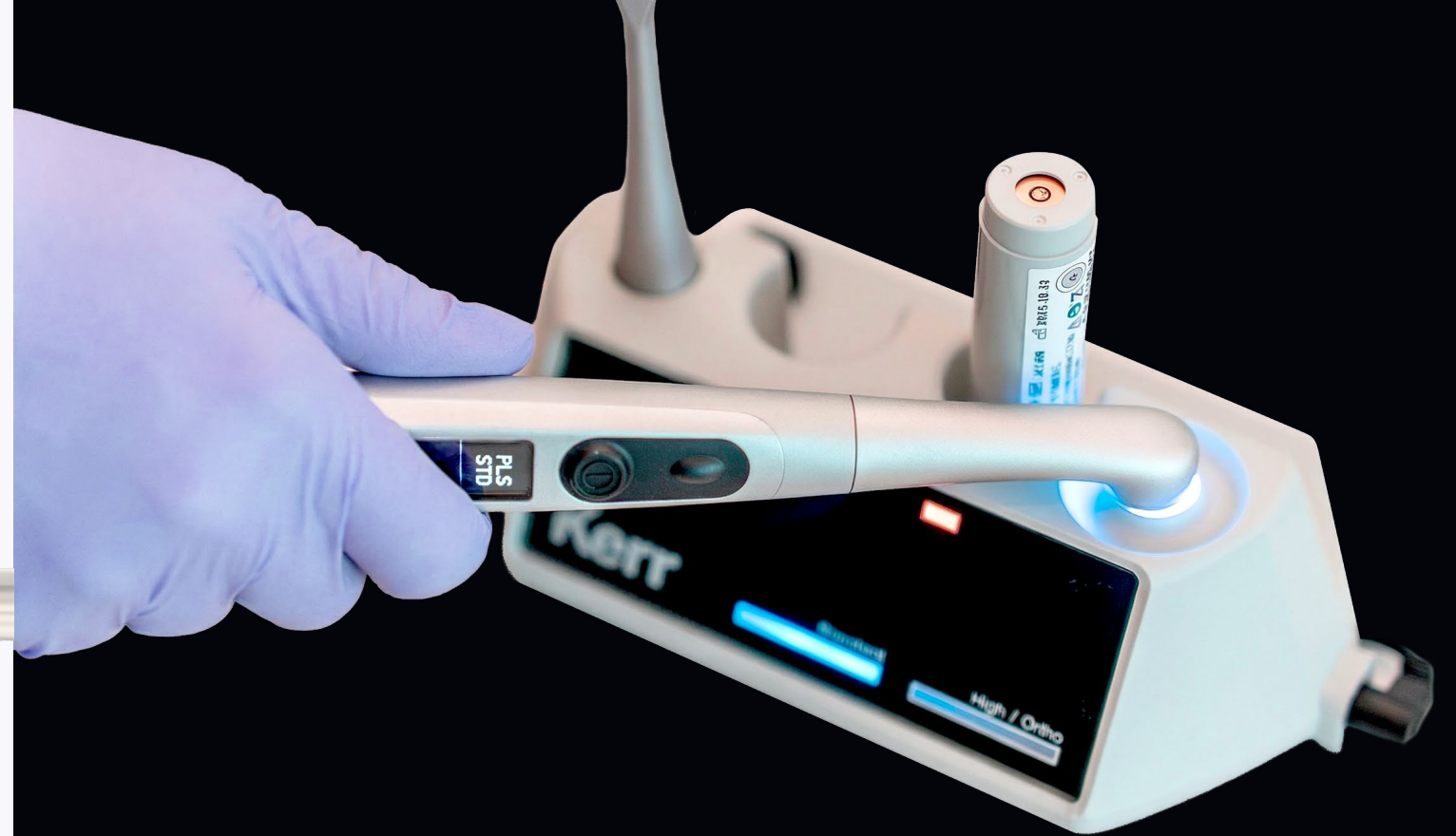
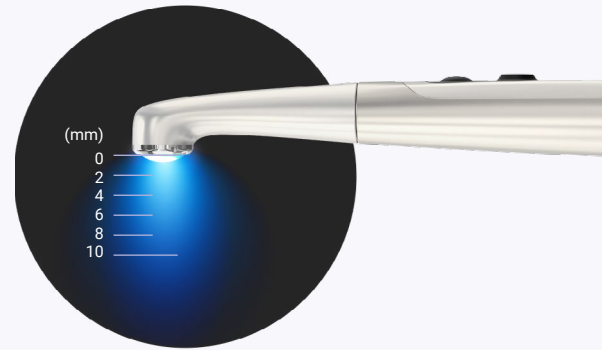
Irradiance at Distance

Demi Pro - Standard Mode (base & pulse average)






Deeper Cures. Powerful Results.

Delivers powerful, collimated light output even at increased distances, ensuring reliable curing performance across various clinical scenarios.



Deeper Cures With Consistency

Demi Pro delivers reliable curing results with consistent, high-intensity light distribution at clinical distances.

			
	Demi Pro	VALO™ X	SmartLite® Pro
Mode	Standard	Std Power (Steady Blue)	Standard
Test Result (0mm) Light tested directly on sensor	1393 mW/cm ²	1137 mW/cm ²	1233 mW/cm ²
Test Result (6mm) Light tested using tooth	1025 mW/cm ²	830 mW/cm ²	748 mW/cm ²

* Data on file

Confidence Starts With Safety

Demi Pro includes a built-in radiometer to easily test its output, giving you peace of mind going into every procedure.

Demi Pro employs the power of Periodic Level Shifting, (PLS), shifting output intensity multiple times throughout the curing cycle to drive continuous high output without generating excess heat.

Demi Pro with PLS 5-Second Exposure



Transillumination Tip*

A powerful tool to aid in intraoral visualization, diagnosis, and treatment planning.



Interproximal caries



Diagnosis of a vertical fracture and occlusal decalcification



Diagnosis of a vertical crack on a tooth undergoing endodontics

Ordering Information

Kit

- 37370 Complete Kit**
1x Curing Tip, 1x Handpiece, 1x Charging Base, 2x Batteries, 3x Light Shields, 1x Tack Cure Lens, 1x AC Adapter

Accessories

- 37371 Curing Tip
37372 Handpiece
37373 Charging Base
37374 Battery
37375 Light Shield Refill
37376 Tack Cure Lens
37377 Transillumination Tip
37378 AC Adapter

Efficient, broad-spectrum curing

Low-profile, 360° rotatable curing head

3 curing modes (Standard, High, and Ortho) with multiple curing times

Lightweight, ergonomic, metal design

Transillumination tip* for diagnosis

*Sold separately

Multiple, long-lasting batteries

Built-in radiometer to inductive charging base



Learn More
kerrdental.com/en-us/demi-pro

*Sold separately

Dr. Saad is a paid consultant of Kerr. The opinions expressed in this article are those of Dr. Saad. Kerr is a medical device manufacturer and does not dispense medical advice. Clinicians should use their own professional judgment in treating their patients.

1-800-KERR-123 | kerrdental.com

©2026 Kerr Corporation. All trademarks are property of Kerr Corporation. All Rights Reserved. Data on File. MKT-25-1336 Rev 1

Kerr™