

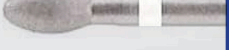


DIAMOND BUR GUIDE



Diamond burs have colored strips at their shank indicating their grit. The greater the coarse, the faster it cuts the tooth. The coarse head is more suitable for cutting objects, and the fine head offers excellent surface finishing and polishing.

UNIVERSAL BUR COLOR CODING FOR STRIPES ON FINISHING AND POLISHING DIAMOND BURS

COLOR and GRIT		PROCEDURES
	BLACK (Supercoarse)	 Crown & Bridge
	GREEN (Coarse)	 Endo Burs
	BLUE (Medium)	 Crown & Bridge Burs
	RED (Fine)	 Composites
	YELLOW (Super Fine)	 Composites (Flow)
	WHITE (Ultra Fine)	 Composites

DIAMOND BURS UNIVERSAL SHAPES GUIDE

BUR TYPE	SURFACES
Flame Burs 	Distal, Mesial, Facial, Buccal and Lingual of Posterior
Football Burs 	Occlusal, Buccal, Lingual anterior

BUR TYPE		SURFACES
Ball (Round) Burs		Occlusal, Lingual anterior
Egg Burs		Occlusal, Buccal, Lingual anterior
Barrel Burs		Occlusal
Inverted Cone Burs		Occlusal
Round End Taper Burs		Distal, Mesial, Buccal, Lingual
Round End Cylinder Burs		Buccal, Lingual
Flat End Cylinder Burs		Buccal, Lingual
Needle Burs		Distal, Mesial, Facial
Wheel Burs		Occlusal, Buccal, Lingual anterior
Depth Cutter Burs		Occlusal

HAVE ADDITIONAL QUESTIONS?

Feel free to contact your Dental City account representative and they'll be happy to help!